



PHILIPS


RESPIRONICS

Amara View

Minimal contact

Full face mask

High marks for **performance**,
no marks on the bridge of
your nose

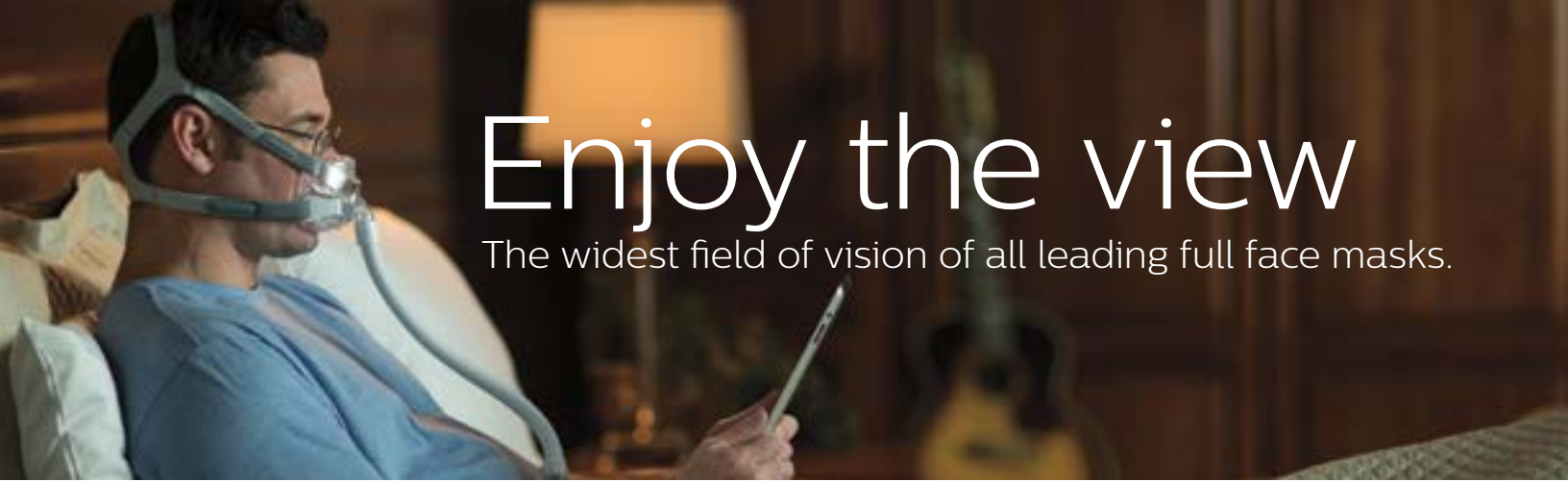


Enjoy the view

The widest field of vision of all leading full face masks.

An idea whose time has come. The Amara View's incredible under-the-nose design makes it the only full face mask that won't cause red marks on the bridge of your nose when compared to the market leaders. Amara View is also the smallest, lightest[†] and has the widest field of vision of all leading full face masks. It's easy to wear glasses, read, watch TV, and use a computer or tablet before falling asleep.*

Features

- Innovative design prevents red marks, discomfort or irritation on the nose bridge*
- The smallest and lightest[†] of all leading full face masks
- Widest field of vision of all leading full face masks*
- Enables wearing glasses, reading, and watching TV
- Easy to use with one modular frame for all cushion sizes
- No bulky cushion and frame in front of your face and eyes

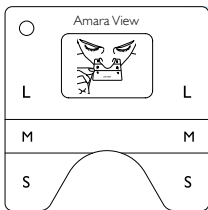

Unlike most full face masks, Amara View's cushion size is determined by measuring the depth of the nose using the simple sizing gauge.



Amara View covers less of a patient's face than all leading full face masks*

With its amazingly innovative design, only Philips Respironics can offer the smallest and lightest choice with the widest field of vision of all leading full face masks.

Get started

To learn more about how Amara View can start making a difference for your patients and your business, contact your Philips Respironics representative for more information.

Item	Part number
Amara View Mask Small	1090602
Amara View Mask Medium	1090603
Amara View Mask Large	1090604

Replacement parts	Part number
Amara View Quick-Release Tube	1090689
Amara View Cushion Small	1090692
Amara View Cushion Medium	1090693
Amara View Cushion Large	1090694
Amara View Headgear Clip 2pk	1090698
Amara View Headgear Large	1090696
Amara View Headgear Standard	1090697

*Internal testing; data on file

† Without quick release tube

Amara View fitting guide

Achieving the right fit with Amara View is easy; just follow the instructions below to get started.
Please read OIS for further fitting instructions.

1



Loosen headgear to a large setting, then grip and twist the headgear clips away from the frame to disconnect.

2



Hold the mask against your face. Place the mask under your nose and over the cushion nostril opening. Push the mask cushion up.

(Note: Do not place the mask over your nose or place your nose in the cushion nostril opening.)

3



Pull the mask over your head.

4



Press the headgear clips into place to connect to mask.

5



Pull back the tabs on the headgear to evenly adjust the top and bottom straps. Pull back the tab to adjust the crown strap.

6



Connect the flexible tubing (included with your CPAP or bi-level device) to the Amara View quick release tube.



©2015 Koninklijke Philips N.V. All rights reserved.
Specifications are subject to change without notice.
Trademarks are the property of Koninklijke Philips N.V.
(Royal Philips) or their respective owners.

Broudy JS 10/12/15 MCI 4106977 PN 1125775

Philips Respironics
1010 Murry Ridge Lane, Murrysville, PA 15668 USA

www.philips.com/respironics